bee2gether vibe™

bee2gether vibe™ Operating Instructions

PLEASE READ CAREFULLY BEFORE OPERATING

Introduction

The bee2gether vibe™ is a vibrating couples ring that is uniquely designed to enhance women's sexual pleasure during intercourse and other types of sexual play. It can also be used alone for self-pleasuring.

When used for intercourse, the bee2gether vibe™ is worn by the man as a “ring” over the penis. Many men like to wear the bee2gether vibe™, for the unique feeling of tightness and engorgement that it provides. The bee2gether vibe™ can contribute to firmer erections, and also to delaying and intensifying orgasmic sensations. **A small amount of water based lubricant is recommended for maximum comfort when putting on, wearing, and taking off the bee2gether vibe™.**

The bee2gether vibe™ should be worn only during use, and for a maximum of 20 minutes at a time.

Components

Start-up and Charging

Your bee2gether vibe™ comes in an attractive carrying case that doubles as both a storage and charging unit. The bee2gether vibe™ will not be fully charged when you receive it. It is recommended that you fully charge the bee2gether vibe™ before initial use. Initial charging may take up to 3 hours.

Firmly place the bee2gether vibe™ with the ring facing upwards and the pleasure ticklers facing forwards into the rear opening of the case. The remote sits facing forward in the front opening of the case, as shown below.

The case has a USB charging port on the back end of the base of the case. Plug the USB cable (provided) into the USB port and plug the other end of the USB cable into the wall power adapter.

Alternatively you can plug the charging base directly into a USB port on your computer or other electronic device with a charging USB port.


When the bee2gether vibe™ is charging, the clear ring on the bottom of the charging case will emit a blue light pulse. The blue light will continue to pulse until the bee2gether vibe™ is fully charged. When fully charged the blue light will glow a constant blue light. Charging time will depend on the battery drainage at the time that it is placed in the case, varying from 5 minutes to up to 3 hours (for a fully drained battery).

The bee2gether vibe™ and remote can be safely stored in the charging case without causing the battery to degrade.


Using the bee2gether vibe™

The bee2gether vibe™ is most commonly worn with the ring around the base of the penis and the stimulator directed upward with the “ticklers” facing forward (i.e. toward partner).

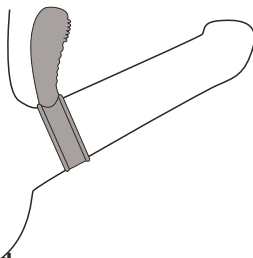For maximum comfort water based lubricant is recommended.

Slide the bee2gether vibe™ over the penis until it sits comfortably against the scrotum. The bee2gether vibe™ glides on most easily over an erect penis, but you can also put it on over a flaccid or semi-erect penis. A small amount of a water-based lubricant can be used to help glide the bee2gether vibe™ over the penis. Once the bee2gether vibe™ is comfortably in place, activate the unit (i.e. switch it to “on” position) by pressing the ON/OFF button on backside of the bee2gether vibe™ (as shown below). The ON/OFF button should be held down for 1.5 seconds to turn the unit on or off.

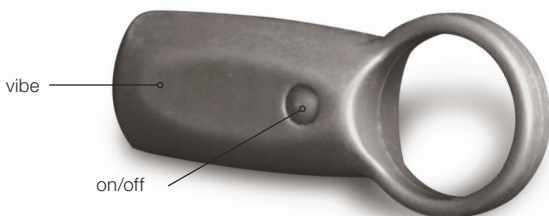Once the bee2gether vibe™ has been activated, the wireless remote can then be used to change the speed and mode as desired (see instructions for wireless remote below).


When using the bee2gether vibe™ for intercourse, the vibe will be ideally touching and stimulating the woman's "pleasure zone" (on and around the clitoris), providing very arousing sensations. Feel free to use the bee2gether vibe™ to experiment with various positions of intercourse and sexual play to find those that are most arousing to you and your partner.

To remove the bee2gether vibe™ after use, turn off the motor (by pressing the ON/OFF button on the backside of the bee2gether vibe™ for 1.5 seconds), grasp the ring and gently glide it off slowly. It may also be helpful to apply a small amount of water-based lubricant to assist with removal. The bee2gether vibe™ can be removed when the penis is in an erect, semi-erect, or flaccid state.

The bee2gether vibe™ has an automatic shut-off feature that turns the unit off after 5 minutes of inactivity in the lowest speed setting. Note that the automatic shut-off is only active when the bee2gether vibe™ is on its lowest speed setting – In all other settings, use the ON/OFF button on the backside of the bee2gether vibe™ to shut the unit down. This automatic shut-off safety feature prevents accidental battery drain if the user puts the bee2gether vibe™ in its lowest speed setting and then forgets to turn it off. Additionally, the bee2gether vibe™ will automatically shut off (from its lowest speed setting) when it is placed in its storage/charging case.

Wireless Remote Control

NOTE: Before using the wireless remote control, you must first activate the bee2gether vibe™ by pressing the **ON/OFF** switch on the backside of the bee2gether vibe™ for 1.5 seconds.

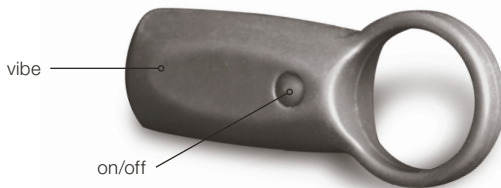There are (3) buttons on the remote control.

The buttons on the remote control function as follows:

- Press the **MODE** button to cycle the bee2gether vibe™ through 4 pre-set vibration sensations
- Press the **+** button to increase the intensity of the bee2gether vibe™ vibration sensations
- Press the **-** button to decrease the intensity of the bee2gether vibe™ vibration sensations

For optimal use, the remote should be held between hip and shoulder height, and pointed towards the vibe.

Cleaning Instructions

Turn the bee2gether vibe™ off prior to cleaning (as described on page 8). The bee2gether vibe™ is waterproof and can be easily cleaned with mild soap and warm water.

To clean the bee2gether vibe™:

1. Turn off the unit
2. Wet the unit by running it under the tap for 5-8 seconds
3. Wash the surface of the unit with a mild soap
4. Rinse under the tap for 5-8 seconds
5. Wipe the unit dry with a clean dry cloth or towel.

Unlike the bee2gether vibe™ unit, the wireless remote is **not** waterproof. **Do not** immerse the remote in water as doing so can damage the remote. To clean the wireless remote, use a slightly damp cloth to wipe it clean. Dry the remote with a clean cloth or towel.

Be sure that both the bee2gether vibe™ and the remote are completely dry before placing them back in the storage/charging case.

Safety and Warnings

Use caution when applying, wearing, and removing the bee2gether vibe™. Be gentle and use slow movements when gliding the ring on and off the penis to avoid unnecessary discomfort (e.g. pulling of pubic hair).

Do not wear/use the bee2gether ring for more than a **MAXIMUM of 20 minutes**. Wearing any stimulator ring for extended periods of time can cause soreness and damage to nerves and tissues. Allow at least one hour between repeated episodes of use.

If you experience pain at anytime while using the device STOP use immediately

- If you notice any swelling, pain, numbness, coldness, bruising, or discoloration of the genitals, remove the stimulator ring **IMMEDIATELY** and consult a physician.
- Never attempt to open any of components of the bee2gether vibe™.
- Do not force the power supply plug into an electrical socket.
- To avoid electric shock, keep the charging case safely away from water and liquids at all times, and use only the bee2gether AC adapter intended for this product.
- Turn “off” after use.
- Do not use the bee2gether vibe™ while charging.
- Only use water soluble lubricants with the bee2gether vibe™.
- Always clean the bee2gether vibe™ after use (follow cleaning directions as described above).
- Do not immerse the hand-held remote in water or other liquids.

- Store in a cool, dry place.
- Never wear the bee2gether vibe™ while sleeping.
- Never use the bee2gether vibe™ while operating a moving vehicle or any machinery.
- Do not use the bee2gether vibe™ on the neck area.
- Do not insert the bee2gether vibe™ into the anus or rectum.
- Seek medical attention immediately if you experience a prolonged erection of more than 3-4 hours.
- Do NOT use the bee2gether vibe™ if:
 - > you have any pain, swelling, tenderness, or bruising on or around the genitals.
 - > you have peripheral vascular disease or any condition that affects normal sensation or blood flow.
 - > you have any condition that interferes with your ability to properly operate and use this device as outlined above.
 - > you have used drugs or substances that could seriously impair your judgment, sensation, or level of alertness.
 - > you have an open sore, infection, or inflammation in the genital area.

Consult your physician prior to use if you have any medical conditions, or if you are uncertain about your ability to safely and/or comfortably use the bee2gether vibe™.

Warranty

The one-year warranty covers product replacement for manufacturing defects only. Visit bee2gethervibe.com for product replacement and warranty information. No other warranties are expressed or implied. If there is a warranty claim, or any other damage, the user's sole remedy is product replacement. No medical claims are implied or warranted by use of this product.

Compliance

bee2gethervibe™ remote: Model 986908.

FCC ID: 2AFYM-0142015, IC: 20660-0142015

bee2gethervibe™ is rated for indoor use from +10C to +45C. This class B digital apparatus complies with Canadian ICES-003. This device complies with part 15 of the FCC Rules and with the industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

bee2gethervibe™ couples vibrator: Model 986908.

FCC ID: 2AFYM-0142015, IC: 20660-0142015

bee2gethervibe™ charging base: Model 986908 Part #9083.

This device complies with part 18 of the FCC Rules and Industry Canada RSS 216.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. For further compliance information, consult bee2gethervibe.com.

Responsible Disposal

The crossed-out wheeled bin on the product indicates that this product should not be treated as household waste. Electrical and electronic equipment contains hazardous substances that can be harmful to the environment and human health if not disposed of properly. To reduce the impact on the environment, we encourage you to use the appropriate take-back systems when the product has reached end of life. Please contact your local or regional waste administration for details.


Contact us for more information on the environmental performance or our products.

Warning

Sold as an adult novelty only; not for medical use. This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities. To avoid electric shock, do not use charger near water. Do not use on swollen or inflamed areas of skin. Discontinue use if discomfort occurs. Do not use on neck area.

bee2gethervibe.com

bee2gether vibe™

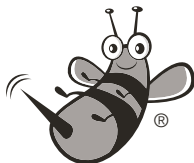THE bee2gether vibe™ IS SOLD AS AN ADULT NOVELTY ONLY. IT IS NOT FOR MEDICAL USE AND HAS NOT BEEN CLINICALLY TESTED. NO MEDICAL CLAIMS ARE WARRANTED OR IMPLIED BY USE OF THIS PRODUCT. THE SUPPLIER BEARS NO RESPONSIBILITY FOR THE PRODUCT SHOULD IT BE MISTREATED, MISUSED, OR USED FOR PURPOSES OTHER THAN AS A PERSONAL MASSAGER.

All materials used to manufacture this product meet RoHS, CE, and/or conform to FDA guidelines.

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